
Instructions for Use

COLD RUSH[®] COMPACT

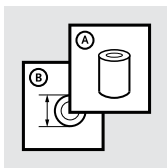3

EN | INSTRUCTIONS FOR USE



5

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11

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18

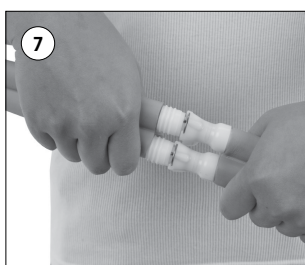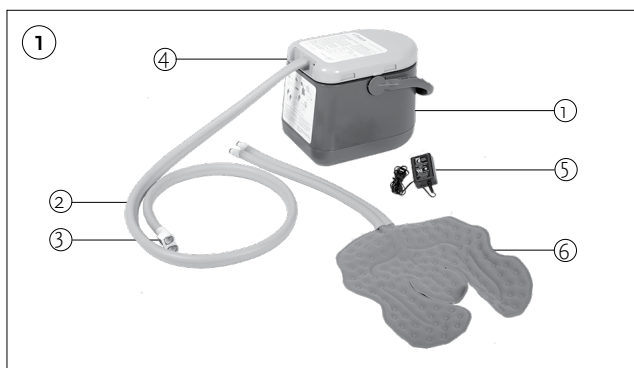EN - Caution: This product has been designed and tested based on single patient usage and is not recommended for multiple patient use. If any problems occur with the use of this product, immediately contact your medical professional.

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ATTENTION - See operating instructions

Meets IEC type BF leakage current requirements

IP21

Protected against solid foreign objects of 12.5mmØ and greater. Protected against vertically falling water drops.

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ENGLISH

INSTRUCTIONS FOR PATIENTS AND MEDICAL PROFESSIONALS

Patients and Medical Professionals should carefully read all of the following instructions and warnings before operating the Cold Rush Compact Cold Therapy System.

INDICATION FOR USE:

- The Cold Rush Compact Cold Therapy System is designed to reduce localized pain and swelling following surgical procedures or trauma.
- The Cold Rush Compact Cold Therapy System is intended for single patient use only.

WARNING: RX ONLY.

This device should only be used with a prescription from a physician that includes the following treatment information: (1) the number of days or weeks that the treatment should last and (2) the length and frequency of product use (and breaks) during treatment. **DO NOT USE THIS DEVICE IF YOU DO NOT HAVE A PRESCRIPTION FROM YOUR PHYSICIAN THAT INCLUDES ALL SUCH TREATMENT INFORMATION.**

A blank treatment protocol to be filled out by your healthcare provider is provided on top of the Cold Rush Compact Cold Therapy Device.

WARNING: POSSIBLE COLD-RELATED INJURIES.

This device reduces the temperature of the skin and tissue. If not used properly and in accordance with the prescribed instructions from your physician, this device can cause serious injury, including frostbite, tissue necrosis and nerve damage.

WARNING: SKIN INSPECTION AND ADVERSE REACTIONS.

You or your healthcare provider must check your skin at least every 1-2 hours during use of this device, regardless of whether you are asleep or awake, for any change in skin condition, including for any increase in pain, burning, blistering, itching, increased swelling, discoloration of the skin, increased redness, welts or any other changes in the appearance of the skin. If you experience any of these reactions, you should immediately discontinue use of your Cold Rush Compact Cold Therapy Device and contact your healthcare provider. Any dressing, casting, bracing or wrapping must not interfere with your ability to check your skin condition. If you are unable to check your skin condition for any reason, do not use the Cold Rush Compact Cold Therapy Device.

WARNING: NO DIRECT CONTACT WITH SKIN; INSULATION BARRIER REQUIRED.

The Cold Rush Pad should never make direct contact with the skin. The Cold Rush pad operating temperature is too cold to be placed directly on the skin. An insulation barrier between the pad and the skin should be used at all times without exception so as to prevent the skin from becoming too cold. Failure to use an insulation barrier between the pad and the skin can lead to a cold-related injury.

WARNING: USE ONLY ÖSSUR COLD RUSH PADS.

Competitive brand cold therapy pads should not be used with the Cold Rush Compact Cold Therapy Device as the temperature may become too cold when using competitive pad and the pads may not otherwise

function properly, causing injury. Only Össur Cold Rush Pads should be used with the Cold Rush Compact Cold Therapy Device.

WARNING: NO MOISTURE IN INSULATION BARRIER.

If the moisture is present in the insulation barrier, the skin may become colder than intended. The insulation barrier should be inspected regularly to check for moisture caused by bleeding, sweating or condensation. If moisture is found on the insulation barrier between the pad and the skin, immediately discontinue use of the Cold Rush until the moisture is removed and a fresh barrier is placed between the pad and the skin.

CONTRAINDICATIONS:

The Cold Rush Compact Cold Therapy Device should not be used on patients with the following conditions:

- History of cold-related injury or adverse reactions to cold therapy;
- Unconsciousness or incoherence due to general anesthesia, sedation or coma;
- Impairment of circulation or potential impairment of healing in the treatment area, including infection;
- Circulatory disease, including Reynaud's disease, Burger's disease, peripheral vascular disease, vasospastic disorders, sickle cell anemia, or hypercoagulable clotting disorders; or
- Polyneuropathy or other nerve damage.

RISK FACTORS:

Physicians should carefully consider the following conditions or factors before prescribing the Cold Rush Compact Cold Therapy Device:

- Sensitivity to cold;
- General impaired circulation;
- Use of medications that may negatively affect peripheral vascular circulation, including beta adrenergic blockers and local epinephrine use (such as local anesthetics);
- Poor nutrition, smoking, tobacco use, excessive use of caffeine or alcohol and any other behavior negatively affecting circulation;
- Desensitization of the treatment area due to local anesthesia or regional nerve blockers;
- Cognitive impairment and use of medications that have a negative effect upon mental capacity or judgment;
- Moisture at the application site due to excessive bleeding, sweating or condensation;
- Diabetes;
- Use of product on hand, wrist, foot or ankle; and
- Young children and elderly.

DEVICE PARTS IDENTIFICATION (Figure 1):

1. Cold Rush Compact Cold Therapy System
2. Hose
3. Couplings
4. DC Connector 12V Power Supply
5. 12V Power Supply
6. Cold Rush Pad (Sold Separately)

APPLICATION/ OPERATION:

1. Ensure cold therapy device is turned off. Do not plug in to wall socket.
2. Fill reservoir with water to the indicated fill line (**Figure 2**).
3. Add chopped or cubed ice to the reservoir up to the indicated fill line

(Figure 3).

4. Close the reservoir with lid and lock with handle. Attach the provided 12V power supply to DC power port, do not plug into wall circuit **(Figure 4).**
5. Apply insulation barrier over the patient's skin in the area to be treated **(Figure 5).**
6. Apply Cold Rush Pad over insulation barrier. The Cold Rush Pad should never be applied directly over the patient's skin **(Figure 6).**
7. Connect Cold Rush pad to hose extending from the device. To ensure proper connection, push couplings together until they snap into place **(Figure 7).**
8. Plug the provided 12v power supply to wall socket (120V US) **(Figure 8).**
9. To turn off the cold therapy device, unplug the power supply unit from wall socket.
10. Always turn the cold therapy device off before disconnecting the Cold Rush Pad or any hoses. To disconnect the hoses, press the release tabs on the coupling and pull the connectors apart.
11. Drain all of the water out of the cold therapy device after each use.
12. Follow directions 2–8 before each use.

HEALTHCARE PROVIDER RESPONSIBILITIES:

- The prescribing physician must determine (1) the number of days or weeks that the treatment should last and (2) the length and frequency of product use (and breaks) during treatment. All such information should be included in the patient's prescription.
- Appropriate training must be provided by the healthcare provider, in the proper application, use and operation and care of the Cold Rush Compact Cold Therapy System.
- The healthcare provider must monitor the patient's use of the device to assure compliance with the prescribed protocol, appropriate use, proper application and operation of the Cold Rush Compact Cold Therapy System, including but not limited to application and maintenance of an insulation barrier between the patient's skin and the Cold Rush Pad.

Treatment Time Period		Frequency and Duration of Use	Inspect Skin Every
Day	Through		
Day			
Day	Through		
Day			
Day	Through		
Day			

ADDITIONAL WARNINGS:

- The Cold Rush Compact Cold Therapy Device should never be left

unattended when plugged in.

- Do not place the tubes or the power cord where your or other's feet can get entangled, causing a fall.
- Never drop or insert any object into any opening or hose.
- Do not operate where aerosol (spray) products are being used or where oxygen is being administered.
- Do not attempt to drag or reposition the Cold Rush Compact Cold Therapy Device by pulling or tugging on the insulated hoses, as this may cause damage to the device and/or cause the hoses to disconnect. Use handle when moving the Cold Rush Compact Cold Therapy Device.
- Do not place the tubes or the power cord by the top of your bed where they could twist around your neck while you are sleeping.
- The pump in Cold Rush Compact Cold Therapy Device is designed to run with water. Running the device without water will cause permanent damage to the pump.
- A wireless communications equipment such as wireless home network devices, mobile phones, cordless telephones and their base stations, walkie-talkies can affect the Cold Rush Compact Cold Therapy Device and should be kept at least a distance 3.3 m away from the Cold Rush Compact Device.

SAFETY STANDARDS:

Cold Rush Compact Cold Therapy System is tested and certified to comply with the IEC60601-1 standard of electrical safety of medical devices and IEC/EN60601-1-2, electromagnetic compatibility for medical electrical devices. IEC 60601-1-11, requirements for medical electrical equipment and medical electrical systems used in the home healthcare environment. The company fulfills the requirements of ISO13485.

ELECTRIC SHOCK HAZARD:

To create the safest operating conditions possible, Össur has designed the Cold Rush Compact Cold Therapy System to operate on a 12 volt power supply, significantly reducing the chance of dangerous electric shock. Avoid use near any water source.

- Only use the power supply provided with the device.
- Unplug device before filling the it with ice and water.
- To avoid danger of electric shock be certain that hands are dry before inserting or removing the power supply from wall socket.
- Avoid getting water or ice on the power adapter.
- Always unplug this product immediately after each use.
- Never use Cold Rush Compact Cold Therapy device while bathing.
- Do not place or store product where it can fall or be pulled into a tub or sink.
- Do not place in or drop into water or other liquid.
- Do not reach for a product that has fallen into water. Unplug immediately.
- Keep cord & device away from heated surfaces.
- Never operate this product if:
 - The device has a damaged electrical cord or plug;
 - The device has been dropped or damaged;
 - The has been dropped into water; or
 - The device is not working properly.
- If any of the above have occurred, immediately discontinue use of the Cold Rush Compact Cold Therapy Device and return the product to

your healthcare provider for examination and repair.

TECHNICAL SPECIFICATION: ELECTRICAL - POWER SUPPLY

- Model no: SA 1312–120100 US
- Input: 120V US
- Input frequency: 60 Hz US
- Output power: 12W
- Output Voltage: 12V
- Output Current: 1A

MAINTENANCE

Do not attempt to repair your Cold Rush Compact Device as this may create a user hazard and will automatically void the warranty.

CLEANING AND PROLONGED STORAGE

- When you have finished all of your cold therapy treatments, prior to prolonged storage, fill the device so the water level is at the ice fill line.
- Add 2 tablespoons of liquid bleach or 3 ounces of 3% hydrogen peroxide to the water and mix.
- Connect hoses to the Cold Rush Pad, plug in the device to start.
- Allow water mixture to circulate for 5 minutes.
- To turn device off, unplug power supply.
- Carefully drain the water mixture out of the device and store in a dry, cool, and dark place leaving the lid slightly open. (Exposure to sunlight and extreme heat may damage hoses and device.)
- When reusing the device after storage, fill the reservoir with clean water and turn the pump on to circulate the water through the hoses and flush the remaining water mixture out of the system (use old pad for this procedure, discard the pad after cleaning is completed).
- Dispose water and refill the reservoir with fresh ice and water combination to begin cold therapy treatment. Follow operating instructions 2–9 before each use.
- Before you use the device, inspect hoses and new pad for any tears or breaks. Contact your healthcare provider to purchase new Cold Rush Pad.
- Clean or replace insulation barrier as recommended by your healthcare practitioner.

ENVIRONMENTAL

- Environmental Operating Temperature: +10°C (+50°F) to 40°C (104°F)
- Operating Humidity: 30% - 75% Relative Humidity
- Atmosphere Pressure: 700 – 1060 h Pa
- Shipping and Storage Temperature: -20°C (-4°F) to 70°C (158°F)
- Shipping and Storage Humidity: 10% - 90% Relative Humidity, non- condensing
- Atmosphere Pressure: 700 – 1060 h Pa

DISPOSAL

Follow local ordinance or regulations for proper disposal of device, accessories, and packaging.

ACCESSORIES

The following accessories are available for your Cold Rush Compact Cold Therapy System at an additional charge. To order, contact your healthcare provider.

- Cold Rush Ankle Pad
 - Cold Rush Hip Pad Left
 - Cold Rush Hip Pad Right
 - Cold Rush Knee Pad
 - Cold Rush Lumbar Pad
 - Cold Rush Rectangular Pad Small
 - Cold Rush Shoulder Pad Small
 - Cold Rush Shoulder Pad Large
 - Cold Rush Universal Pad
- Insulation barrier used for testing thermal properties of the system is a single layer of 3-ply Seal-In Sock, Össur product numbers (K-531310, K-531313, K-531316, K-531320).

LIMITED WARRANTY

Össur warrants, to the original purchaser, that this product is free from defects in materials and workmanship for 6 months from initial purchase. This warranty does not apply if this product: 1 was not purchased from Össur or its authorized dealer 2 has been altered in any way; or 3 has not been used in accordance with Össur's Instructions For Use (IFU). This Limited Warranty does not cover damage due to accidents, neglect, misuse, or operation beyond capacity, parts damaged by improper installation, substitution of parts not approved by Össur, or any other alteration or repair by others that, in Össur's judgment, materially or adversely affects the product or part. This product expected service life is 6 months from the date of first use.

For complete warranty guidelines please visit www.ossur.com

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